


JAGA Systems

Portable Wireless Electrophysiology

multichannel · compact · affordable

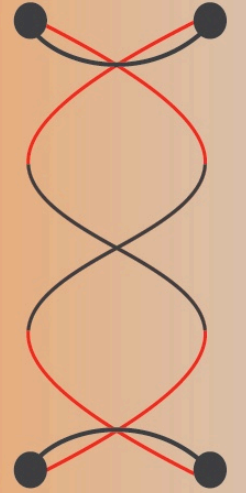


Jinga-hi
Sunnyvale, CA


How to use JAGA16

info@jinga-hi.com <> www.jinga-hi.com <> Sunnyvale, CA


1. What is included in JAGA16

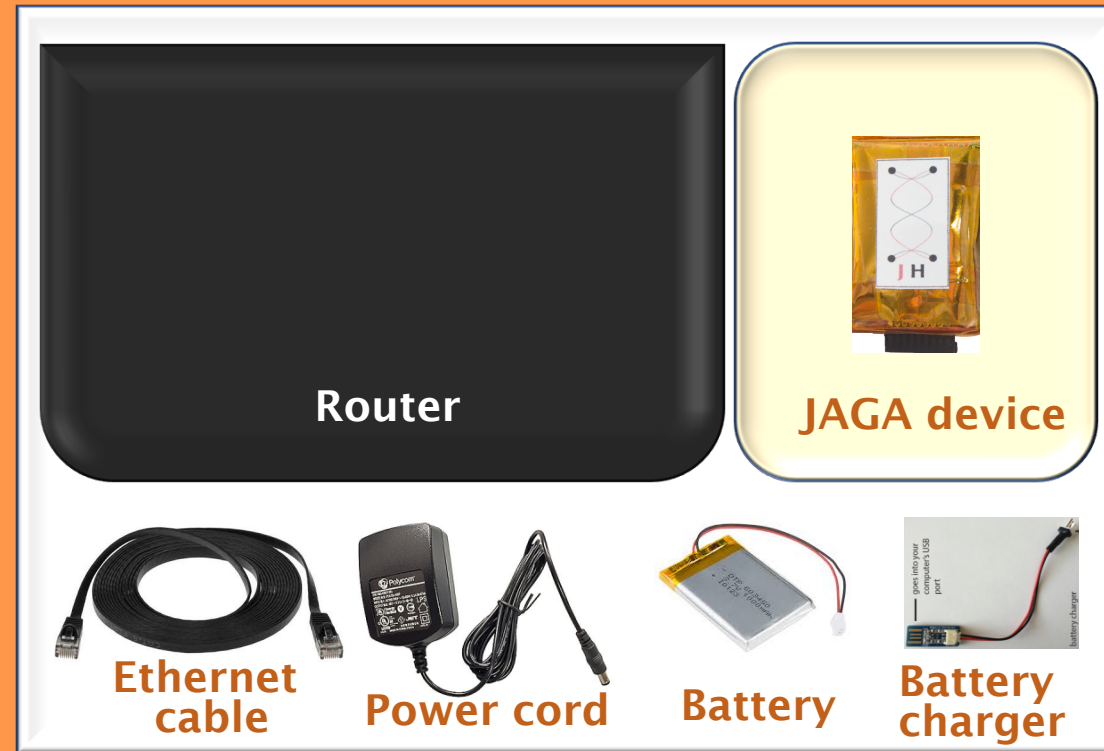

2. Power the router

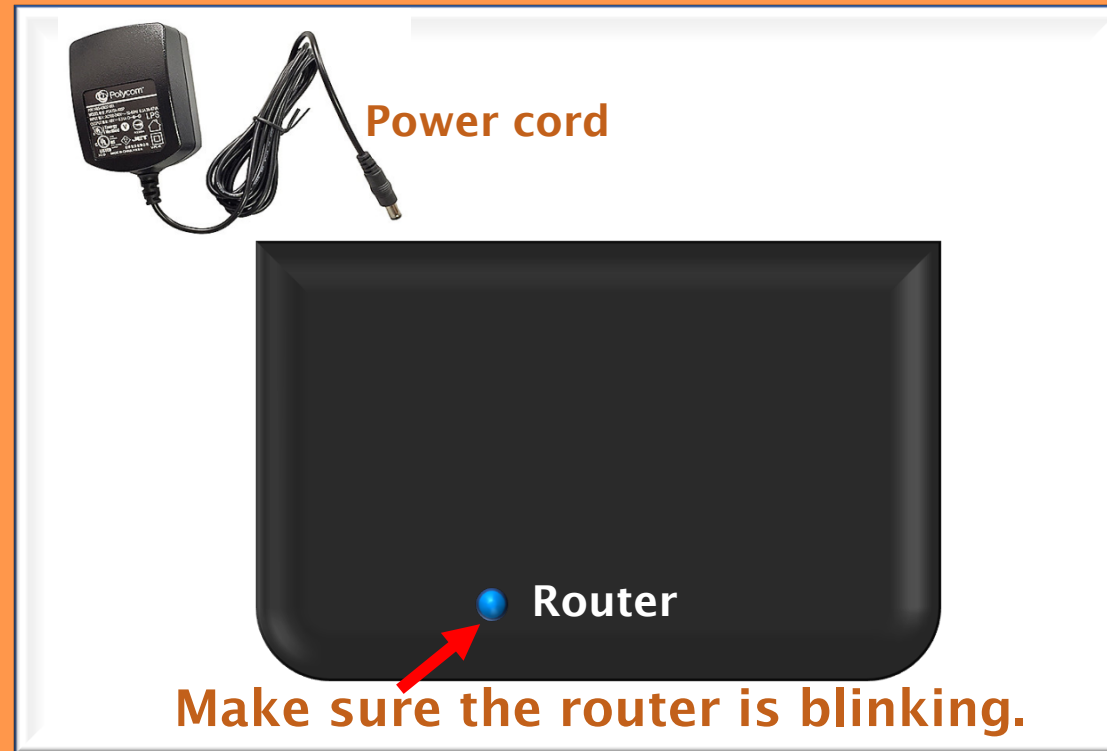


3. Find the jaganet

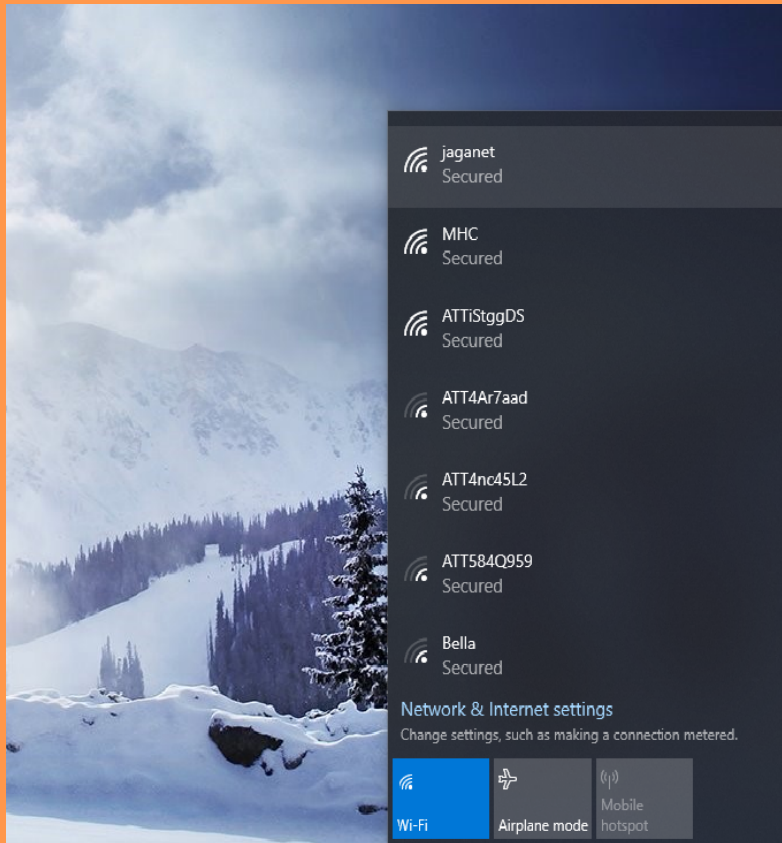

← **jaganet:
network for JAGA
on your computer**

4. Power JAGA

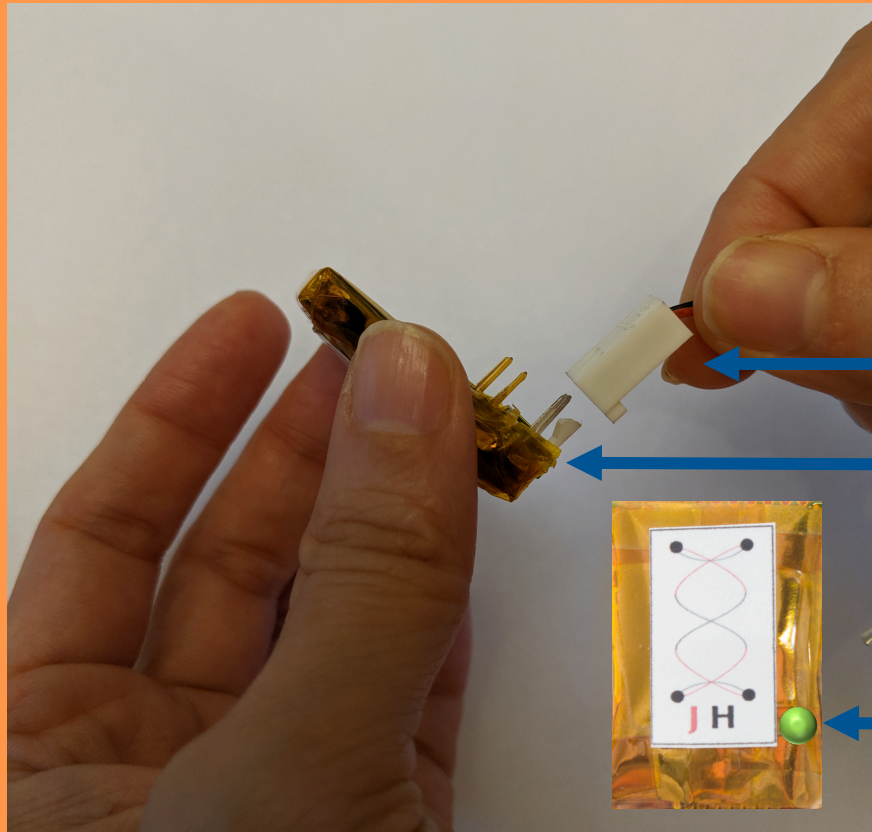


This is battery.

This is JAGA16.


**Green light
should turn on.**

4. Power JAGA

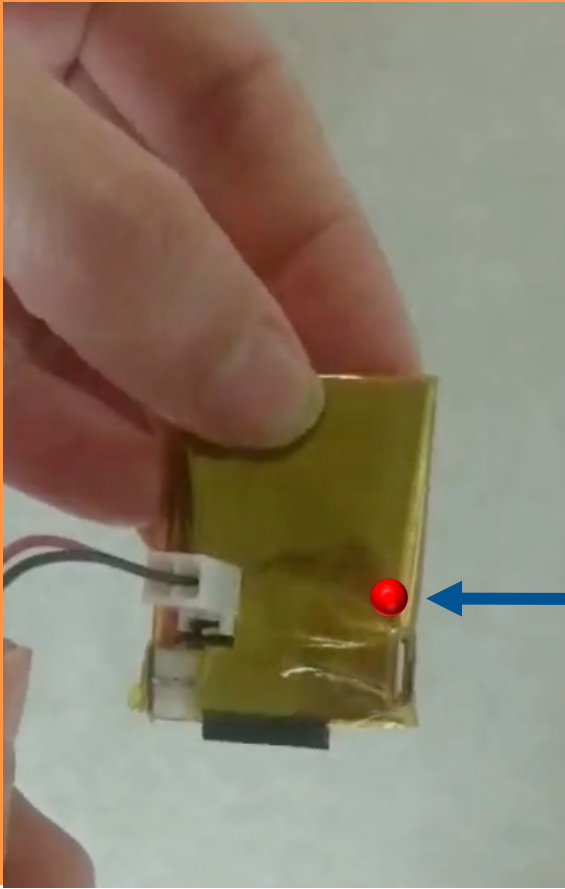


**Red light on the other side
should flash for a few seconds.**

5. Start the software

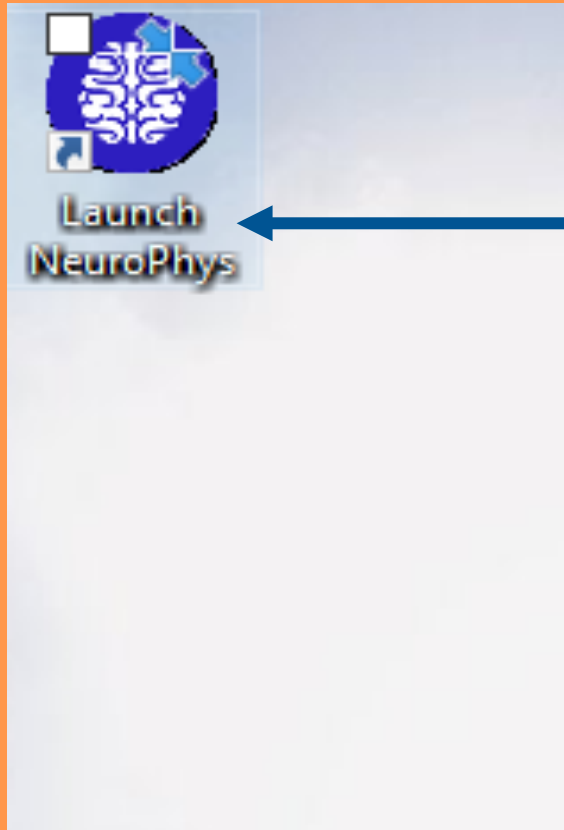


**Install the software first.
Then Launch NeuroPhys**

5. Start the software

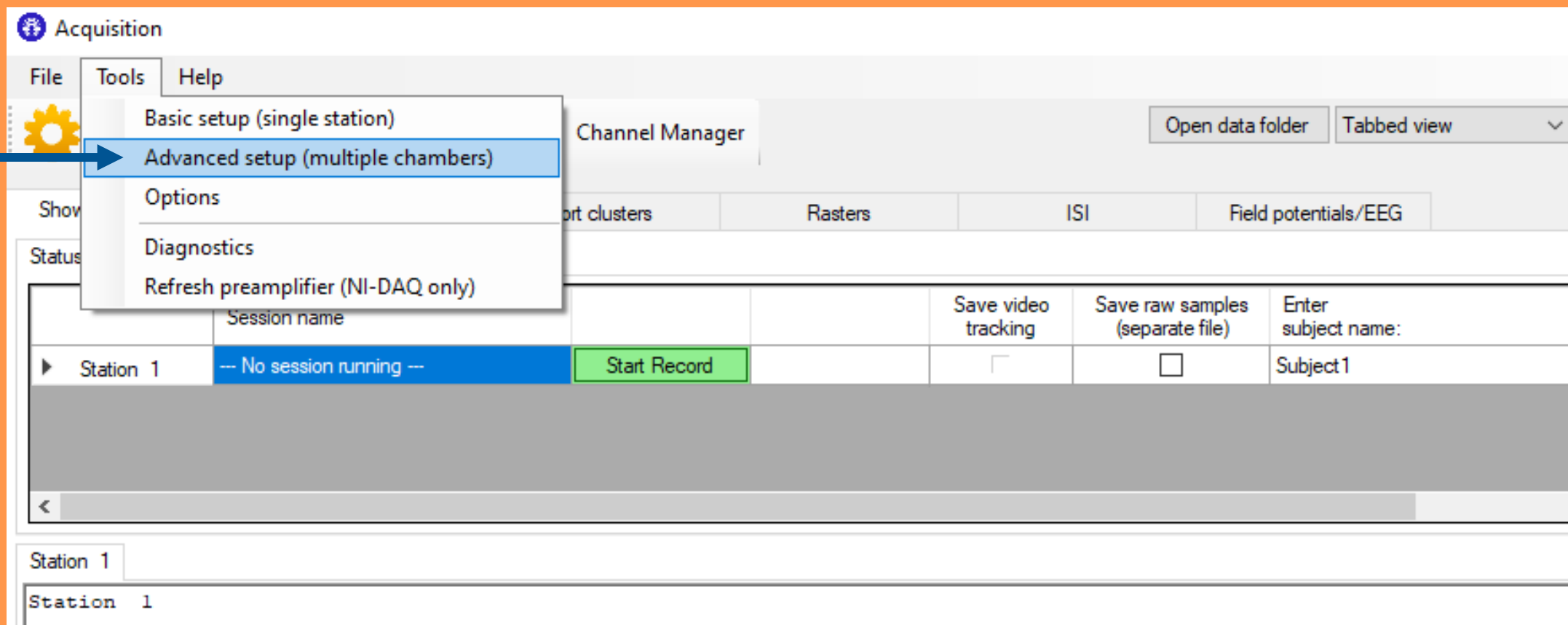

5. Start the software

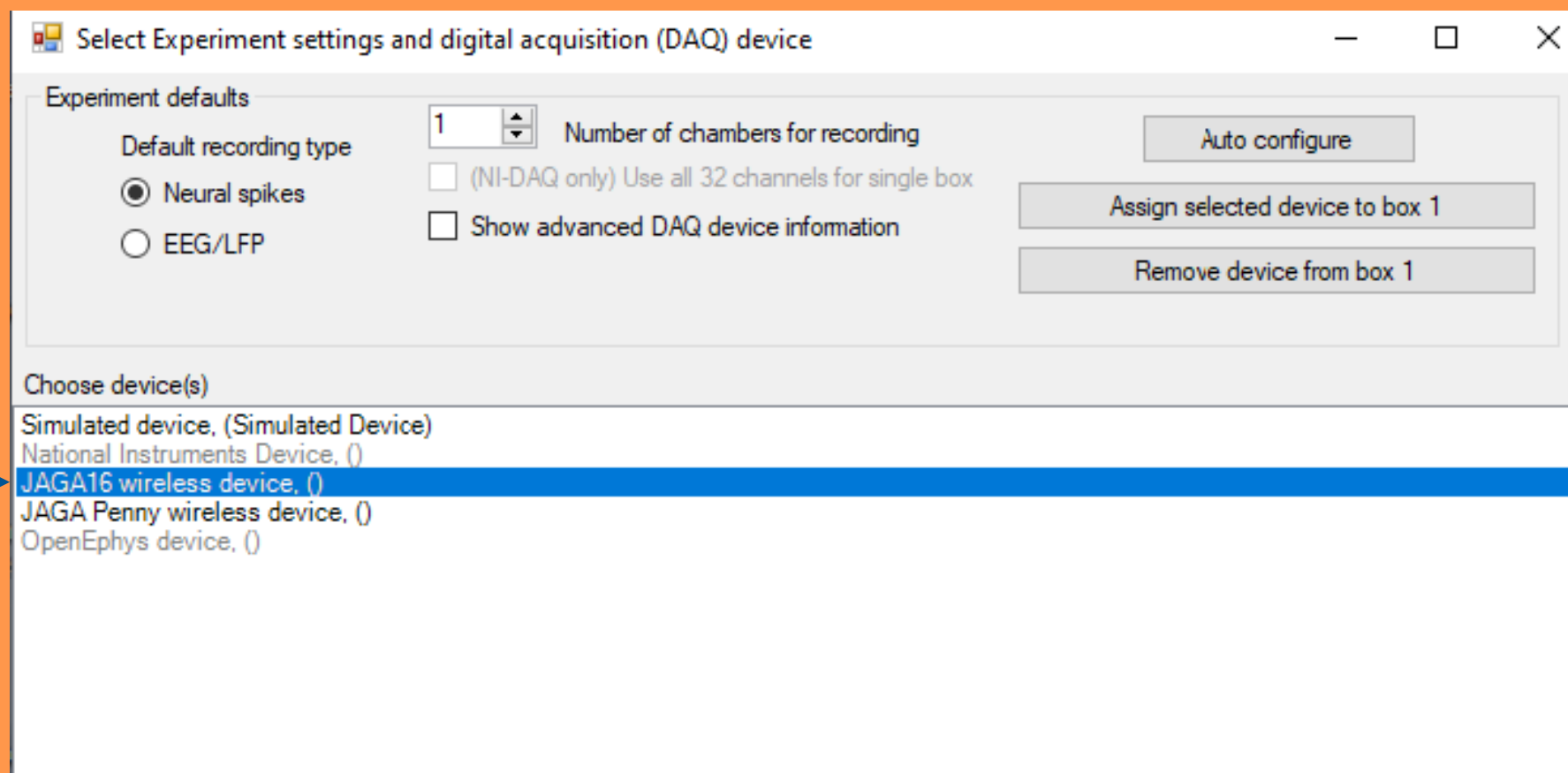


6. Display Data streams wirelessly into your laptop.



6. How to charge battery



**Connect the battery charger
to your computer's USB port.**

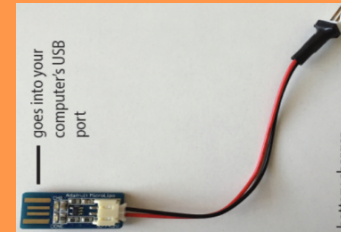

**Connect the battery
to the battery charger.**

Enjoy your experiments!



Thank you for your scientific contribution to the world.

info@jinga-hi.com <> www.jinga-hi.com <> Sunnyvale, CA